
Chiro Vision Upright Radiographic System



A Floor/Wall or Floor/Ceiling mount tubestand with upright wallstand to effortlessly accomplish all standing exams. The stainless steel bearings and electromagnetic locks enable ease of positioning of vertical and horizontal movements. The Chiro Vision is an affordable solution for film based or digital imaging.



**VISION
IMAGING**
PARTNERS, INC.

ChiroVision Upright Radiographic System

Standard Features:

Tubestand

- Mounts floor to wall or floor to ceiling
- Tube rotation enables precision alignment of the tube for non-angled, angled and upright views.
- Convenient display gives indication of tube rotation.
- 6' (183 cm) tracks allow for 40" (100 cm) and 72" (180 cm) SID film distances.
- Stainless steel bearings assure smooth vertical and horizontal tubestand movement.
- Electro-magnetic locks provide positive positioning of vertical and horizontal tubestand positions.

Vertical Wallstand



- Fully counterbalanced for effortless movement.
- 17" X 17" (43 cm X 43 cm) grid cabinet
- Low absorption front panel with patient centering line.
- 103 Line, 8:1 fine line grid for superior film quality.
- Manual cassette tray accepts 14" X 17" (35 cm X 43 cm) and smaller film sizes, both inch and metric.
- Electro-magnetic lock maintains cabinet height.
- Extensive vertical travel of 61.3" (155 cm) allows for views from cervical to weight bearing.

Options

- Choice of high frequency and single phase generators to meet specific needs and budget.
- 17" X 17" (43 cm X 43 cm) bucky for wallstand
- 10:1 or 12:1 ratio grid for wallstand
- 178 line, 10:1 ratio CR grid.
- Size sensing cassette tray for use with PBL collimation.
- Ion chamber to accept AEC capability of generator.
- 14" X 36" (35 cm X 91 cm) grid cabinet or reciprocating bucky for wallstand (Model 3547)

Manufacturer's Statement

- Complies with all CDRH performance standards in effect at date of manufacture.
- Specifications are subject to change without notice.
- Made in U.S.A.
- An ISO13485 & ISO9001 Certified Facility.
- ETL Listed.

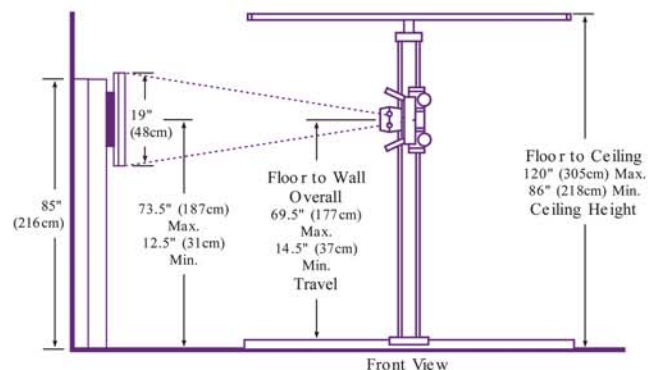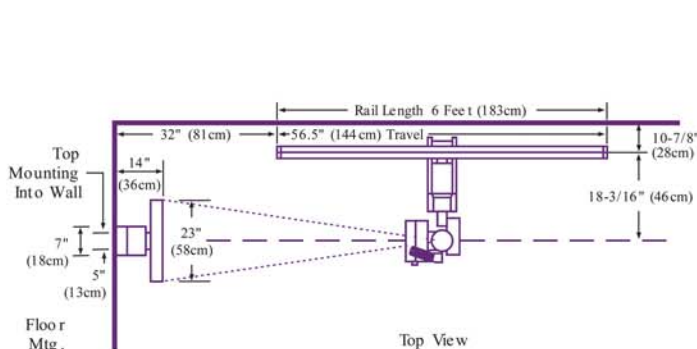

P.O. Box 1425, Hunt Valley, MD 21030 Ph. 856.424.8661 Fax: 856.751.8678